

HEM-7156
& HEM-7156-A

OMRON



Effortless Clinical Accuracy from the Comforts of Home.

3 key benefits of the IntelliWrap™ cuff

360°
Accuracy
Intelli
Wrap™



ACCURATE

360° accuracy with IntelliWrap™ Cuff, every time.



COMFORTABLE

Pre-formed cylindrical cuff fits snugly around upper arm, for all arm sizes between 22-42cm.



INTUITIVE


Easily wearable with just one hand.

*Based on the results of an independent quantitative survey conducted by MIMS in Jul-Sep 2019 on total 600 doctors (including cardiologists, physicians and specialists) in India, Indonesia, Philippines, Malaysia and Vietnam who have recommended BPM products to their patients.

All for Healthcare


HEM-7156 & HEM-7156-A

DID YOU KNOW?

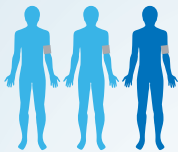

1 in 3 people measure their blood pressure incorrectly¹, while using their blood pressure monitor at home.

Common Mistakes When Wearing Blood Pressure Monitor Cuffs

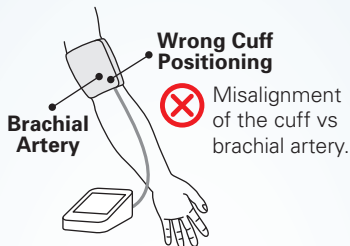

Are you perhaps one of them?

When it comes to blood pressure measurements, cuff positioning greatly contributes to the accuracy of readings. OMRON HEM-7156 blood pressure monitor eliminates this concern. Designed with innovative IntelliWrap™ Cuff, the device delivers consistent and accurate measurements regardless of the cuff's positioning around the upper arm.



Care for your Heart Health

Take control of your heart health and achieve clinical grade accuracy effortlessly at home²⁻³ with OMRON's HEM-7156.


3 Key Benefits of the IntelliWrap™ Cuff




Reduce Errors

Regular cuffs may cause inaccurate results when positioned incorrectly. IntelliWrap™ 360° accuracy eliminates that concern, ensuring consistent and precise accuracy no matter how you wear the cuff.

REGULAR CUFF

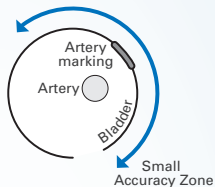


INTELLI WRAP™ CUFF

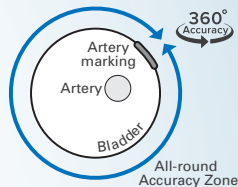


A Fit That Feels Just Right

The IntelliWrap™ is cylindrically pre-formed with tapering ends, allowing it to fit snugly around your upper arm. This secure fit enables not only accurate but also comfortable measurements every time. The wide cuff-size range of 22-42cm makes it suitable for most Asian arm sizes.

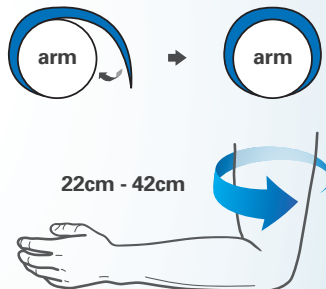

Easy To Wear With One Hand

The IntelliWrap™ is easily wearable with just one hand, perfect for seniors and those who require self-monitoring at home.

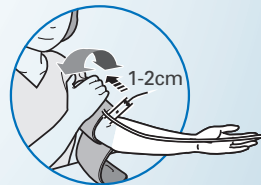

Why Choose OMRON?

**CLINICALLY
VALIDATED**

All Products are Clinically Validated

We ensure 100% OMRON Blood Pressure Monitors are clinically validated for accuracy and meet the highest international standards of precision and reliability.



Intellisense® Technology for Comfort & Accuracy

Intellisense®, available in all OMRON Blood Pressure Monitors, can automatically personalise the cuff inflation just for you to ensure less pain and discomfort during measurement.



Trusted Brand from Japan with Over 45-Years Experience on Home Blood Pressure Monitors

Our blood pressure monitors are available in more than 110 countries and trusted worldwide by millions of patients and healthcare professionals. We offer a wide range of products to suit your needs and keep innovating advanced technology and solution in heart healthcare field.

Specifications for HEM-7156 / HEM-7156-A

Features	IntelliWrap™ Cuff	Body Movement Detection	Pulse
	Intellisense	Hypertension Indicator	Cuff Wrapping Guide
	60 Memory	Average of Last 3 Readings	Irregular Heartbeat Detection
Measurement Method	Oscillometric Method		
Blood Pressure Measurement Range	SYS: 60 to 260 mm/Hg DIA: 40 to 215 mm/Hg		
Pulse Measurement Range	40 to 180 beats/min		
Accuracy	Pressure: ± 3 mm/Hg Pulse: $\pm 5\%$ of reading		
Battery Life	Capacity of new alkaline batteries is approx. 1,000 measurements		
Power Source	4 "AA" alkaline batteries 1.5V or AC adapter (optional, INPUT AC100-240V 50-60Hz, 0.12 - 0.065A)		
Console Weight	Monitor: approx. 340g without batteries • Arm cuff: approx. 163g		
Outer Dimensions	Approx. 105mm (W) x 85mm (H) x 152mm (L)		
Arm Cuff	Fits arm circumference 22 - 42cm		
Package Content (HEM-7156)	Monitor, IntelliWrap™ Cuff, storage case, 4 AA batteries, instruction manual (1) and (2)		
Package Content (HEM-7156-A)	Monitor, IntelliWrap™ Cuff, AC adaptor, storage case, 4 AA batteries, instruction manual (1) and (2)		

References:

1. Perleberg Consumer Research 2014, UK/Italy.
2. Bilo G et al. Impact of cu positioning of blood pressure measurement accuracy. May a specially designed cu make a difference? Hypertens Res 2017. DOI:10.1038/hr.2016.184.
3. IP Protocol, ESH, 2010.

