



# SKB041-3

PATIENT  
TRANSPORTATION  
TROLLEY




**SAIKANG** Medical

Healthier  
your life

# SKB041-3

## FEATURES AND HIGHLIGHTS

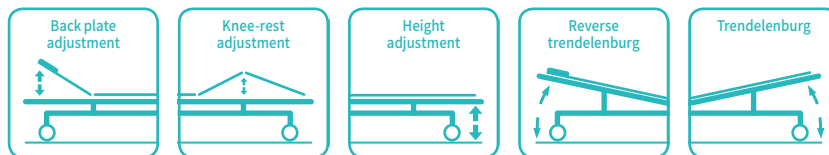

### Technical parameters

▪ External size (LxW):	2100x850 mm
▪ Table height adjustment:	630-930 mm
▪ Back plate adjustment:	0-75° (±10°)
▪ Knee-rest adjustment:	0-45°
▪ Reverse Trendelenburg:	0-15°
▪ Trendelenburg:	0-15°



### Technical configuration

<input checked="" type="checkbox"/> Usa hydraulic pump	2 pcs
<input checked="" type="checkbox"/> Crank system	1 set
<input checked="" type="checkbox"/> Back rest air spring	2 pcs
<input checked="" type="checkbox"/> 8" German castors, central locking	4 pcs
<input checked="" type="checkbox"/> Emergency push handle	1 set
<input checked="" type="checkbox"/> Directional wheel	1 pc
<input checked="" type="checkbox"/> Folding side rail	1 set
<input checked="" type="checkbox"/> Φ74mm wall bumper wheel	4 pcs
<input checked="" type="checkbox"/> Iv pole previsions	4 pcs
<input checked="" type="checkbox"/> Mattress	1 pc
<input checked="" type="checkbox"/> IV pole	1 set
<input checked="" type="checkbox"/> Central locking break pedal	2 sets
<input checked="" type="checkbox"/> Directional wheel release pedal	2 pcs
<input checked="" type="checkbox"/> X-ray translucent bed platform	1 set



**Metal material**

Baosteel from Fortune Global 500  
Panasonic Robotic ensure 360° full  
smooth wedding.

**Plastic**

Environmental plastic can be  
sterilized in 100°C, tensile upto  
30MPa.

**Steels | Painting**

11 Process epoxy painting, ASTM  
testing anti-bacterial, paint thickness  
0.12 mm, brightness 60°, paint can  
resist 50 kg impact.

**Brake**

Break spring keep qualified after  
10,000times uninteruped break.

**Hydraulic pumps**


Height adjustment, trendelenburg  
and revised tendelenburg through  
two USA hydraulic pumps.

**Center wheel**

Can be engaged to facilitate direction  
control, enabling users to navigate  
the stretcher easily in hectic and  
dynamic environments.

**Rising**

Hydraulic treadle for rise.

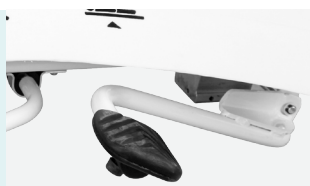**Falling**

① Leg section lows down to make  
head higher.  
② Back section lows down to make  
leg section higher.  
③ Step on the middle of the pedal,  
make the bed base fall down steadily.

**Castors**

Four double side casters of 8 inch  
diameter, simultaneously lockable by  
the step pedal at both sides.  
TPR tire no worn out after running  
30KM, save anti-winding hard shell,  
united forming without bolts.

**Mattress**

Waterproof PU mattress, convenient  
for transfer patient.

**Side rail**

Foldable aluminium alloy guardrail

**Bumper wheel**


Bumper wheel detachable, max.  
protect the bed and patient during  
transportation

**Drainage Hooks**

Mobile plastic drainage hooks.

**IV pole**

Stainless Steel Column with plastic  
hook.

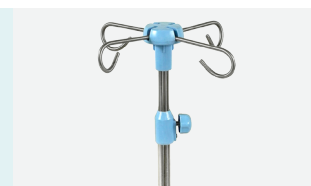



## CONSTANT INNOVATION



Design



Research



Manufacturing



Maintenance



Service



## WORLDWIDE BRANCHES

### South America Branch

**Saikang Medical Technology, Inc**  
W3400 Corrientes, Argentina  
TEL: +55 (11) 2758-0438  
Rua Universal, Hollywood,  
São Bernardo do Campo,  
Brazil  
TEL: +55 (11) 2758-0438

### East Africa Branch

**Saikang Medical Technology Ltd**  
Mezzanine Floor, Taj Tower,  
upper hill, nairobi, Kenya  
TEL: +254714084806

### South East Asia Branch

**Saikang Medical Technology Ltd**  
Sala Thammassop Subdistrict  
Thawi  
Watthana District  
Bangkok Thailand  
TEL: +66 642741046

### UK Branch

**Uk Healthward International Limited**  
DEPT 111,  
196 High Road Wood Green,  
London, N22 8HH  
United Kingdom  
TEL: +31 (0) 202111106

### Russia Branch

**LLC 'SAIKANG'**  
454004, Russian Federation,  
Chelyabinsk,  
Universitetskaya Naberezhnaya,  
Russia  
TEL: +7 (919)300-00-06

### Mid East Branch

**Saikang Medical Supplies (Mid East) LLC**  
Musaffah  
P.O Box 53200, Abu Dhabi,  
UAE  
Tel: +971 585710652

# SAIKANG / Medical

Jiangsu Saikang Medical Equipment Co., Ltd.  
STOCK CODE: 870098

📍 No. 35, Lehong Road · Zhangjiagang (Jiangsu) · China

☎ +86 18021231901 · 🕒 24/7

✉ export@saikangmedical.com

🌐 www.saikangmedical.com



Company video



General catalogue